REV 1

Chip Pak **High Volume Pressure** Sensor for Disposable, **Backside Pressure Applications**

Motorola has developed a low cost, high volume, miniature pressure sensor package which is ideal as a sub-module component or a disposable unit. The unique concept of the Chip Pak allows great flexibility in system design while allowing an economic solution for the designer. This new chip carrier package uses Motorola's unique sensor die with its patented, piezoresistive implant technology, along with the added feature of on-chip, thin-film temperature compensation and calibration.

Features:

Low Cost

- Patented piezoresistive strain gauge implant, temperature compensation and calibration all integrated on a single, monolithic sensor die.
- Pressure Range Available: 0-300 mmHg
- Polysulfone (Mindell S-1000) Case Material (Medical, Class VI Approved)

side of the Chip Pak (pressure is applied to the backside of the device). Front side die and wire protection must be provided in the customer's housing. Use caution when handling the devices during all processes.

Pin Number, Style 1

2

S+

Motorola is offering the Chip Pak option package. Application-specific parts will have an "SPX" prefix, followed by a four digit number, unique to the specific customer. Devices will be shipped in a tape and reel packaging.

NOTE: The die and wire bonds are exposed on the front

MAXIMUM RATINGS			
Rating	Symbol	Value	Unit
Overpressure (Backside)	P _{max}	125	PSI
Burst Pressure	Pburst	1000	kPa
Supply Voltage	V _{Smax}	10	Vdc
Storage Temperature	T _{stg}	85	°C
Operating Temperature	ТА	+15 to +40	°C

1

٧s

Motorola's MPX2300D Pressure Sensors. Motorola's MPX2300D pressure sensor has been designed for medical usage by combining the performance of Motorola's shear stress pressure sensor design and the use of biomedically approved materials. Materials with a proven history in medical situations have been chosen to provide a sensor that can be used with confidence in applications, such as invasive blood pressure monitoring. It can be sterilized using ethylene oxide. The portions of the pressure sensor that are required to be biomedically approved are the rigid housing and the gel coating.

The rigid housing is molded from a white, medical grade polysulfone that has passed extensive biological testing including: tissue culture test, rabbit implant, hemolysis, intracutaneous test in rabbits, and system toxicity, USP.

A silicone dielectric gel that has been used extensively in implants covers the silicon piezoresistive sensing element. The gel is a nontoxic, nonallergenic polymer system which passes pyrogen testing, as well as meeting all USP XX Biological Testing Class VI requirements. The properties of the gel allow it to transmit pressure uniformly to the diaphragm surface, while isolating the internal electrical connections from the corrosive effects of fluids, such as saline solution. The gel provides electrical isolation sufficient to withstand defibrillation testing, as specified in the proposed Association for the Advancement of Medical Instrumentation (AAMI) Standard for blood pressure transducers. A biomedically approved opaque filler in the gel prevents bright operating room




3

S–

4

Ground

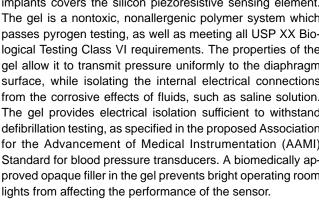





MPX2300D



MPX2300D

OPERATING CHARACTERISTICS ($V_S = 6 \text{ Vdc}, T_A = 25^{\circ}\text{C}$ unless otherwise noted)

Characteristics	Symbol	Min	Тур	Max	Unit
Pressure Range	POP	0	—	300	mmHg
Supply Voltage ⁽⁸⁾	٧ _S	—	6.0	10	Vdc
Supply Current	۱ _۵	—	1.0	—	mAdc
Zero Pressure Offset	Voff	-0.75	—	0.75	mV
Sensitivity	—	4.95	5.0	5.05	μV/V/mmHg
Full Scale Span ⁽¹⁾	VFSS	2.976	3.006	3.036	mV
Linearity + Hysteresis	—	- 2.0	_	2.0	
Temperature Effect on Sensitivity	TCS	-0.1	_	+0.1	%/°C
Temperature Effect on Full Scale Span	TCVFSS	-0.1	_	+0.1	%/°C
Temperature Effect on Offset ⁽⁴⁾	TCV _{off}	-9.0	—	+9.0	μV/°C
Input Impedance	Z _{in}	1800	—	4500	Ω
Output Impedance	Z _{out}	270	—	330	Ω
R _{CAL} (150 kΩ) ⁽⁹⁾	RCAL	97	100	103	mmHg
Response Time ⁽⁵⁾ (10% to 90%)	t _R	—	1.0	_	ms
Temperature Error Band	_	0	—	85	°C
Stability ⁽⁶⁾	—	—	±0.5	—	%VFSS

MECHANICAL CHARACTERISTICS

Characteristics	Symbol	Min	Тур	Мах	Unit
Weight (Case 423)	_	_	170	_	mg
Warm–Up	_	_	15	_	Sec

NOTES:

1. Measured at 6.0 Vdc excitation for 100 mmHg pressure differential. V_{FSS} and FSS are like terms representing the algebraic difference between full scale output and zero pressure offset.

2. Maximum deviation from end-point straight line fit at 0 and 300 mmHg.

3. Slope of end-point straight line fit to full scale span at 0°C and +85°C relative to +25°C.

4. Slope of end-point straight line fit to zero pressure offset at 0°C and +85°C relative to +25°C.

5. For a 0 to 300 mmHg pressure step change.

6. Stability is defined as the maximum difference in output at any pressure within POP and temperature within +10°C to +85°C after:


a. 1000 temperature cycles, -40°C to +125°C.
b. 1.5 million pressure cycles, 0 to 300 mmHg.

7. Operating characteristics based on positive pressure differential relative to the vacuum side (gauge/differential).

8. Recommended voltage supply: 6 V ± 0.2 V, regulated. Sensor output is ratiometric to the voltage supply. Supply voltages above +10 V may induce additional error due to device self-heating.


9. Offset measurement with respect to the measured sensitivity when a 150k ohm resistor is connected to VS and S+ output.

PACKAGE DIMENSIONS

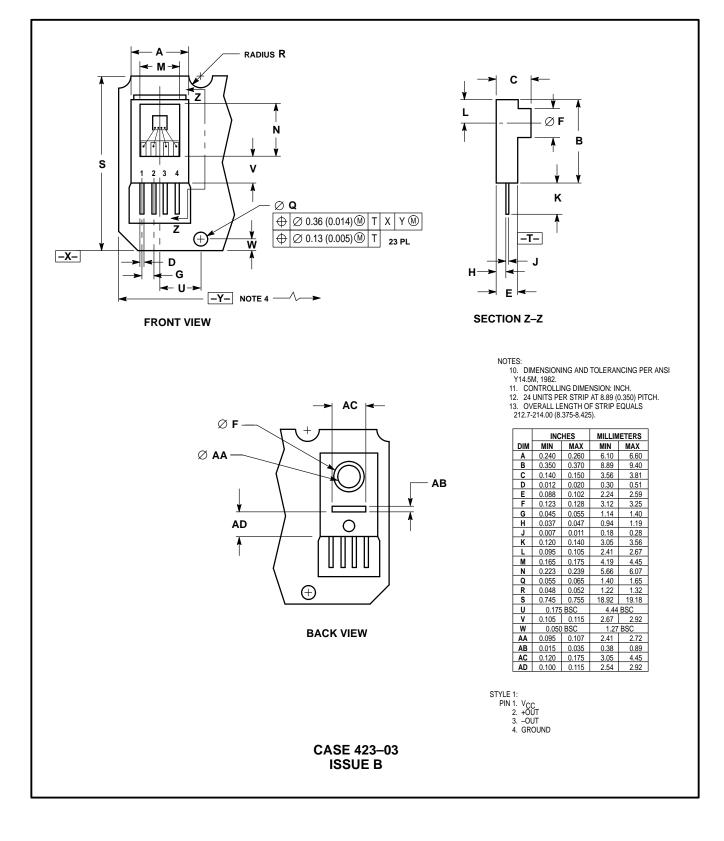

MPX2300D

PACKAGE DIMENSIONS — CONTINUED

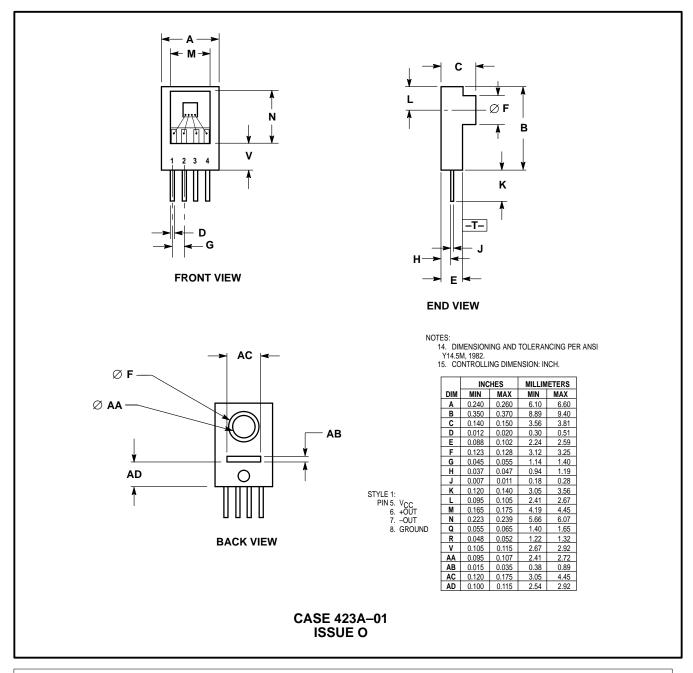

Motorola reserves the right to make changes without further notice to any products herein. Motorola makes no warranty, representation or guarantee regarding the suitability of its products for any particular purpose, nor does Motorola assume any liability arising out of the application or use of any product or circuit, and specifically disclaims any and all liability, including without limitation consequential or incidental damages. "Typical" parameters can and do vary in different applications. All operating parameters, including "Typicals" must be validated for each customer application by customer's technical experts. Motorola does not convey any license under its patent rights nor the rights of others. Motorola products are not designed, intended, or authorized for use as components in systems intended for surgical implant into the body, or other applications intended to support or sustain life, or for any other application in which the failure of the Motorola product could create a situation where personal injury or death may occur. Should Buyer purchase or use Motorola products for any such unintended or unauthorized application, Buyer shall indemnify and hold Motorola and its officers, employees, subsidiaries, affiliates, and distributors harmless against all claims, costs, damages, and expenses, and reasonable attorney fees arising out of, directly or indirectly, any claim of personal injury or death Motorola and its officers. Motorola was negligent regarding the design or manufacture of the part. Motorola and its a negrificant such unintended or unauthorized use, even if such claim alleges that Motorola was negligent regarding the design or manufacture of the part. Motorola and its a negligent regarding the design or manufacture of the part. Motorola and its a negligent regarding the design or manufacture of the part.

Literature Distribution Centers:

USA: Motorola Literature Distribution; P.O. Box 20912; Phoenix, Arizona 85036. EUROPE: Motorola Ltd.; European Literature Centre; 88 Tanners Drive, Blakelands, Milton Keynes, MK14 5BP, England. JAPAN: Nippon Motorola Ltd.; 4–32–1, Nishi–Gotanda, Shinagawa–ku, Tokyo 141, Japan. ASIA PACIFIC: Motorola Semiconductors H.K. Ltd.; Silicon Harbour Center, No. 2 Dai King Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong.



 \Diamond

